

NOSTRESS CHAIR

Manual or electric



The NoStress chair is very comfortable, multi-positional, and facilitates transfers thanks to its removable armrests.

The NoStress can be used both at home and in the communities because of its easy-to-clean cover. It consists of removable, zippered duvet modules that can be washed at 30°.

This recliner has **removable and/or retractable armrests** to facilitate transfers.


The aesthetically pleasing NoStress is an **exclusive evolutionary product** on the market.

You can upgrade your manual NoStress to an electric version at any time.

Features

- **Optimal hygiene:** replaceable, removable and machine washable duvet modules at 30°.
- **Composition of the duvet:** filled with siliconized hollow fibre for ultra-comfort.
- **Practical:** no power outlet required, it can be installed wherever you want, at home or in the community.
- **3 positions:** sit, relax (TV), rest.
- **Safety:** M2 fire-resistant PVC/PU coating.
- **Colours:** 3 coatings available
PVC/PU: Anthracite, Raspberry, Pistachio - **Water-repellent velvet:** Brown or Pearl Grey -
Mottled fabric: Mottled grey
- **Seat width:** 53 cm
- **Height floor/seat:** 48 cm
- **Upgradeable chair:** The manual NoStress chair can easily be upgraded to an electric NoStress
- **Options**
 - Castor kit: consisting of 4 x 5 cm castors (manual model only)
 - Extension kit: consisting of 4 x 8 cm feet (manual model only)
 - Electric kit: for upgrading manual NoStress to electric tilt NoStress
 - HD (High Density) foam
 - Battery Kit
- **Maximum user weight:** 130 kg
- **Warranty:**
 - Chair: 2 years
 - Actuator: 5 years

NOSTRESS remote control electric version

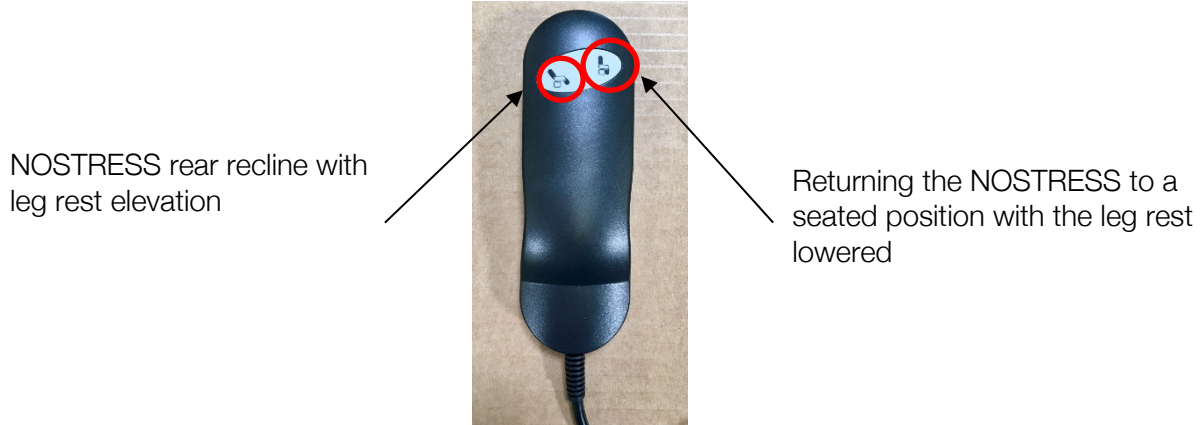

Help with maintenance








The duvet of your chair is removable and can be machine washed at 30°. We advise you to place your duvet in a washing net to preserve its original appearance.

The expected life of the chair is influenced by regular use, storage, maintenance, and cleaning.

To keep your INNOV'SA chair in good condition, you should maintain it regularly as follows:

- Check that screws, levers and adjustment knobs are secure.
- Check the chair for cracks, damage to structural parts, etc. and clean it. Replace the cover if necessary.

Technical specifications

Seat width	Seat depth	Overall width	Seat height	Overall height	Weight of the chair	Maximum user weight
						
54 cm	55 cm	72 cm	48 cm	112 cm	Manual: 43 kg Electric: 47 kg	130 kg
<p>Non-fire standards PVC/PU coating: M2 - 1021-1 - 1021-2 and biocompatible</p> <p>Non-flammable standards for velvet and fabric coatings: 1021-1 - 1021-2 and biocompatible</p>						
<p>Dimensional and weight tolerances ± 15mm / 1.5 kg</p>						