

EU300

TRAÇION-t

4 Section Traction Table (Electric)

EU300 Traction is the ultimate Traction Table for all Traction Needs: Lumbar + Cervical Traction



- 1 Year on Gas Struts, Castors, Upholstery
- 2 Years on Electronic Components
- 3 Years on Frame & Pillar

Prussian Blue

Licorice Black

Dim Grey



2 Motors



Hand Control



Foot Control



2 / 3 Sides Elevating Bar

Maximum User Weight
225 kgs / 496 lbs

Product Features

ALL FOUR CORNERS OF THE EU300 ARE SUPPORTED BY 3 INCH SWIVEL CASTORS, ONE AT EACH CORNER.

- 1st Motorized Actuator for Hi/Lo Adj.
- 2nd Motorized Actuator for Centre Arching Mechanism
- Twin “Oversized” Pillars for Added Stability
- Individual “Leveling Foot” on All 4 Corners
- Retractable Castors at Both Ends
- Twin Gas Struts for Solid Head Section Support
- Maximum Load Capacity of 225 kgs / 496 lbs



Special Features

The EU300 Taction is equipped with ALL NEW “Centre Arching Mechanism” for Motorized LUMBAR ELEVATION. The facilitates any pulling traction force on the lumbar area with a slight upward lifting support (up to 26° RISE) especially when the patient is lying in the supine position.



Dual TSS Switches

For easy accessibility & control, the “Centre Arching/Lifting” is controlled by 2x TSS Switches on both sides just underneath the mid-section of table top.



Adjustable Recoil Spring

The “Adjustable Recoil Spring” helps the 3rd and 4th section of table top gradually returns to its original position as the applied traction force decreases. We implemented a handle at the end of the table (near the 4th section) for recoil spring intensity adjustment. With this feature, therapists can adjust the intensity anytime during the therapy to make the decompression smooth and effective, depending on the patient’s weight.

Ultimate Traction Solution

EU300 Traction + ET800 Traction Unit + CT850 Cervical Traction Add-On Device

the combo are the ultimate solution package for any physiotherapists who are required to perform **Traction and Decompression Therapy** to treat patients suffering from pain and limited function due to lower back and neck conditions.



ET800 Traction Unit (Add-on Device)

Standard Package Contents

- ET800 Traction Unit / Head (Main Traction Unit)
- Power Cord (Plug Type: EU / US / UK / AU / SA)
- Knob Screw for Mounting on EU300 (4 pcs)
- Fuse / Ext. Cord / Manuals / Patient Safety Switch



Versatility

The ET800 Traction Unit can be placed facing the physiotherapist when the table is at the corner of a room. Or it can be placed facing the therapist when operating at the rear-end or when additional room space is available.

CT850 Cervical Traction (Add-on Device)

CT850 Materials & Constructions

With the all-new CT850 Cervical Traction Add-On Device, the neck wedges, head pad and baseboard are all molded from carefully selected high-quality composite polymer. The complete slider is made of light-weight extruded aluminum. Especially the “Block Joints + Pulleys”, they are CNC milled and machined to perfection from a whole billey aluminum.

CT850 Adjustment


The adjustment of traction angle is simply done by the electric or hydraulic feature of the traction table. The neck wedges can be widened or narrowed depending on patient neck width from either side.

CT850 Assembly

All components and parts of the CT850 are hand assembled with precision. Hence from the therapist point of view, installation and dismantling of the CT850 on/off the traction table is simple, swift and effortless. Just pull out the aluminum PIN that is secured by an O-ring and insert the slider's stainless steel PIN to the over-hanger then insert the alu. PIN back one again. Mount the over-hanger on the Traction Unit and secure only ONE knob screw. Then pull the traction cord over the 1st pulley then underneath the 2nd one. Walla, the CT850 is ready to provide the patient with a comfortable cervical traction therapy.



Note: The ET800 Traction Unit(Head), CT850 Cervical Traction Add-on Device, Thoracic Traction Belt and Lumbar Traction Belt shown in the picture herein are **OPTIONAL** and **SOLD SEPARATELY**.

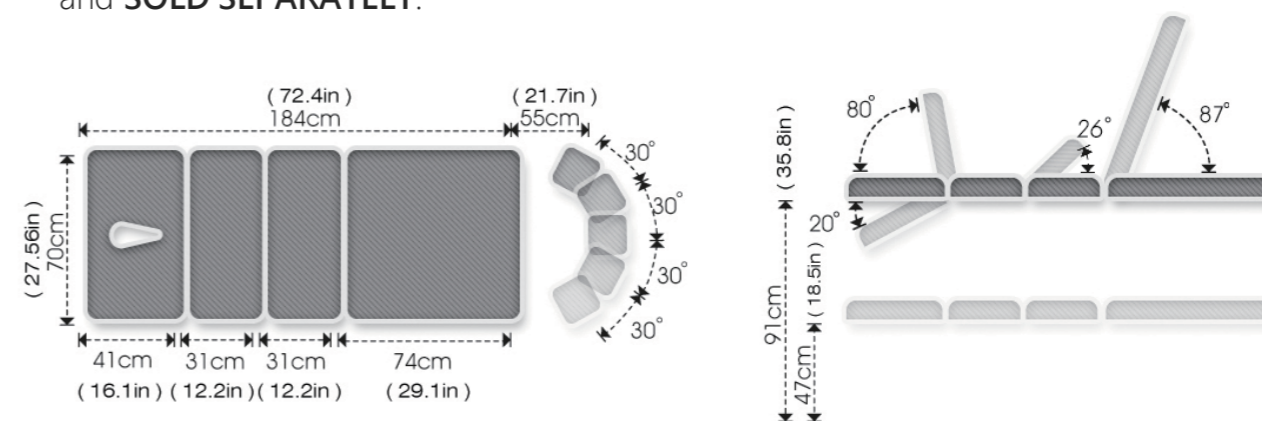

ET800

Traction Unit

TRAÇION-U



• 1 Year Full Manufacturer Warranty

EVERTRAC[®]

Maximum Pulling Force
91 kgs / 200 lbs

Standard Package Contents

- ET800 Traction Unit / Head (Main Traction Unit)
- Power Cord (Plug Type: EU / US / UK / AU / SA)
- Knob Screw for Mounting on EU300 (4 pcs)
- Fuse / Ext. Cord / Manuals / Patient Safety Switch

Product Features

Traction-u ET800 is the result of thousands of feedbacks from tens of thousands of physiotherapists over the past decade. The Retro look of Traç Pod is also one of the many features that physiotherapists requested to retain for good. Instead of a touch screen LCD panel and working thru complicated windows behind windows, physiotherapists recommended the easy-to-adjust “HOLD” and “REST” knobs as they have been proven to be one of the fastest and effective way of adjustment of traction force during treatment. The “HOLD” and “REST TIME” can be adjusted in real-time as well as the FORCE. Once set, press START and the new setting will take effect in the next coming cycle. The unit is equipped with a “Patient Switch” for Full Control of an Emergency Situation, all at the hands of a patient. The Traçion-u ET800 works best on Traçion EU300 however it also can be fitted on tables from other brands.

For Cervical Traction Treatment, CT850 is available as an optional add-on part. To enable this function, a switch behind the unit can be located to switch between (1) Lumbar Traction and (2) Cervical Traction. Weight reading of kgs and lbs can also be switched via a switch behind the unit.

With the ET800, CT850 and EU300, you have the ultimate traction system in your home, clinic, hospital or rehab center.

Versatility

The ET800 Traction Unit can be placed facing the physiotherapist when the table is at the corner of a room. Or it can be placed facing the therapist when operating at the rear-end or when additional room space is available.

Controls Behind

